

Pedi Fit

Tourniquet Cuffs

INVENTED FOR CHILDREN





Pedi Fit

Tourniquet Cuffs



INNOVATIONS FOR CHILDREN

Children are special, and they require special tourniquets for surgery.

Delfi's **Pedi Fit Pediatric Cuffs** with matching limb protection sleeves were invented specifically for children to meet their unique requirements.

Reusable and Sterile Single Use cuffs are available in four cuff sizes to fit the limbs of most children, from neonates to adolescents.

Ideal for use with the Delfi **PTS ii Portable Tourniquet System** and compatible with most automatic tourniquet instruments.

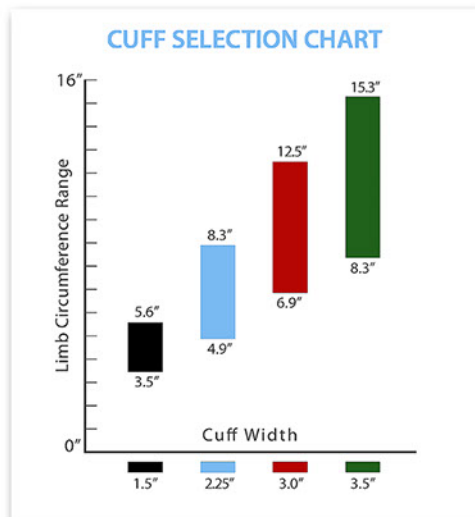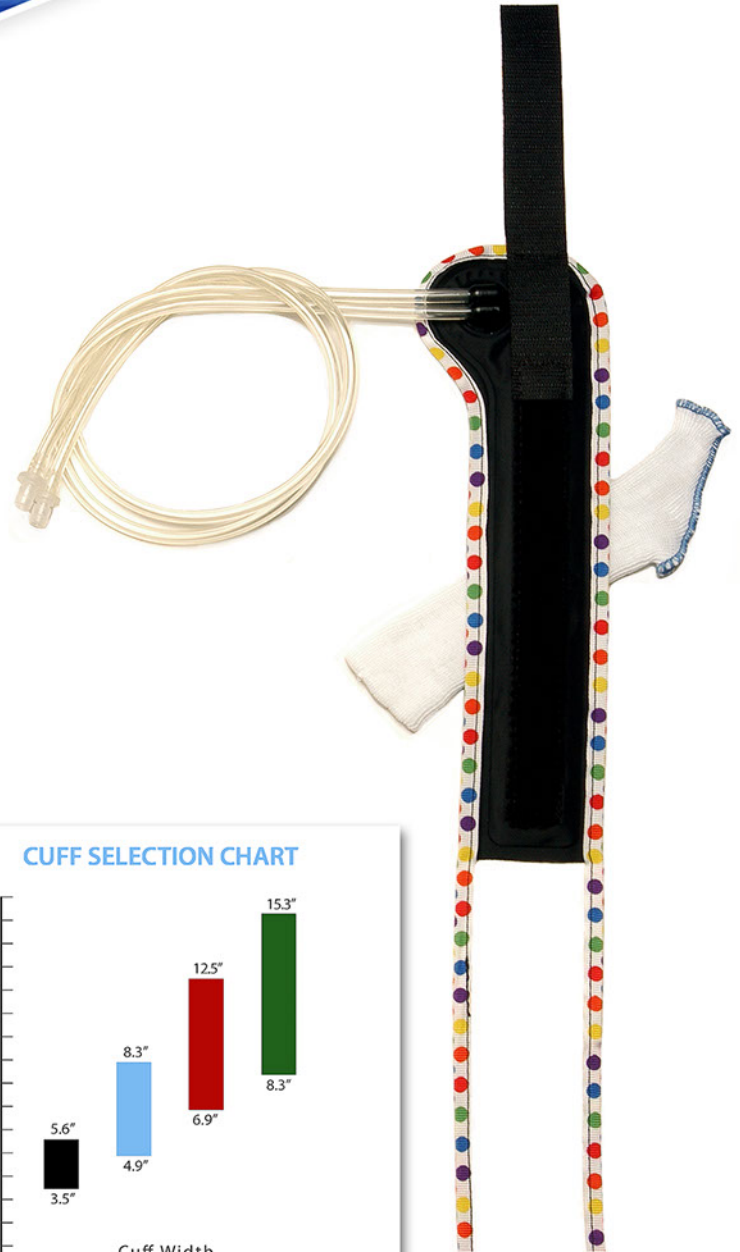



PEDI FIT TOURNIQUETS

HELPING IMPROVE PEDIATRIC SURGERY

* **SOFT & FLEXIBLE**

Design of both cuffs and matching limb protection sleeves allow close conformance around pediatric limbs for the application of tourniquet pressure

* **P-SHAPED CUFF DESIGN**

An offset port allows the cuff to conform to small circumference limbs and guides tubing away from the surgical site, which reduces obstruction and makes draping of the limb and cuff easier

* **TWO INDEPENDENT CUFF FASTENERS**

Help make proper cuff application easier, and provide a secondary means of securing the cuff around the limb

* **OPTIMIZED CUFF WIDTH**

The width of each cuff is designed to reduce the need for high tourniquet pressures while providing the maximum surgical exposure

* **PATENTED POSITIVE LOCKING CONNECTORS**

Fully compatible with Delfi's **PTSii Portable Tourniquet Systems** and A.T.S.* Automatic Tourniquet Systems

MATCHING LIMB PROTECTION SLEEVES

DESIGNED TO HELP PROTECT A CHILD'S LIMB FROM POSSIBLE PINCHING AND BRUISING

* **PROVEN PERFORMANCE**

Delfi Limb Protection Sleeves have been proven to be more effective than other common forms of limb protection in reducing wrinkling and pinching of skin

* **CUFF MATCHED**

Each color coded sleeve is designed to match a specific **PediFit** size and provide optimum limb protection without affecting cuff performance and stability

* **CONVENIENCE**

Easy to apply matching limb protection sleeves are included with all reusable and sterile single-use **Pedi Fit** Tourniquets





Pedi Fit

Tourniquet Cuffs

Available in four sizes from neonates to adolescents.

Delfi Medical Innovations is a world leader in advancing tourniquet technology and safety. Members of Delfi's team have been exclusively involved in tourniquet development and evaluation for more than 30 years. Delfi's patented tourniquet instruments and tourniquet cuffs have been designed to help surgeons and operating-room nurses provide safe and effective care for their patients.

ORDERING INFORMATION

STERILE SINGLE-USE (BOX OF 10)

SINGLE PORT	DUAL PORT	SIZE
9-9700-001	9-9800-001	1.50" W X 7.25" L
9-9700-002	9-9800-002	2.25" W X 10.5" L
9-9700-003	9-9800-003	3.00" W X 15.0" L
9-9700-004	9-9800-004	3.50" W X 17.5" L

Each Sterile Single Use PediFit Cuff is supplied with a matching limb protection sleeve

REUSABLE (1 EACH)

SINGLE PORT	DUAL PORT	SIZE
9-7700-001	9-7800-001	1.50" W X 7.25" L
9-7700-002	9-7800-002	2.25" W X 10.5" L
9-7700-003	9-7800-003	3.00" W X 15.0" L
9-7700-004	9-7800-004	3.50" W X 17.5" L
9-7700-005	9-7800-005	SET OF 4 (EACH SIZE)

Each Reusable PediFit Cuff is supplied with 5 non-sterile matching limb protection sleeves

MATCHING LIMB PROTECTION SLEEVES (10 PER BOX)

9-7900-001	black trim	1.50 sleeve
9-7900-002	blue trim	2.25 sleeve
9-7900-003	red trim	3.00 sleeve
9-7900-004	green trim	3.50 sleeve

CONTACT INFORMATION

P: 800.933.3022 +604.742.0600

F: +604.742.3800

E: info@delfimedical.com

Suite 106 - 1099 West 8th Avenue,
Vancouver, B.C., Canada, V6H 1C3

www.delfimedical.com

U.S. Patent No. 5,741,295; 5,649,954; 5,578,055; 5,454,831 and foreign patents and patent applications

LATEX FREE

D-4900-003 Rev 2 © 2012

 **Delfi** medical
innovations inc.