

Easi-Fit



Variable Contour and
Cylindrical Tourniquet Cuffs

**Unique Easi-Fit design and
patented features provide:**

- * Personalized limb-shape matching
- * Easy, consistent cuff application
- * Superior limb protection matched to each Easi-Fit cuff

Extended ports allow pneumatic connections to be made away from the patient's limb with positive locking connectors.

Offset port design takes the guesswork out of port location, allows the ports to be positioned in the desired location prior to securing the cuff around the limb.




Contour cuff shape in combination with a unique variable fastener allows a superior fit along both proximal and distal edges of the cuff.

Variable fastener allows the cuff to conform to a wide range of limb shapes for uniform transfer of cuff pressure across the entire cuff width to help lower limb occlusion pressure.



Application handle makes cuff application more consistent, opens to make cuff removal faster and helps maintain the cuff in a stable position on the limb. No annoying tie ribbons.

Secondary fastener for enhanced stability and ease of application. (Underside of cuff near variable fastener.)

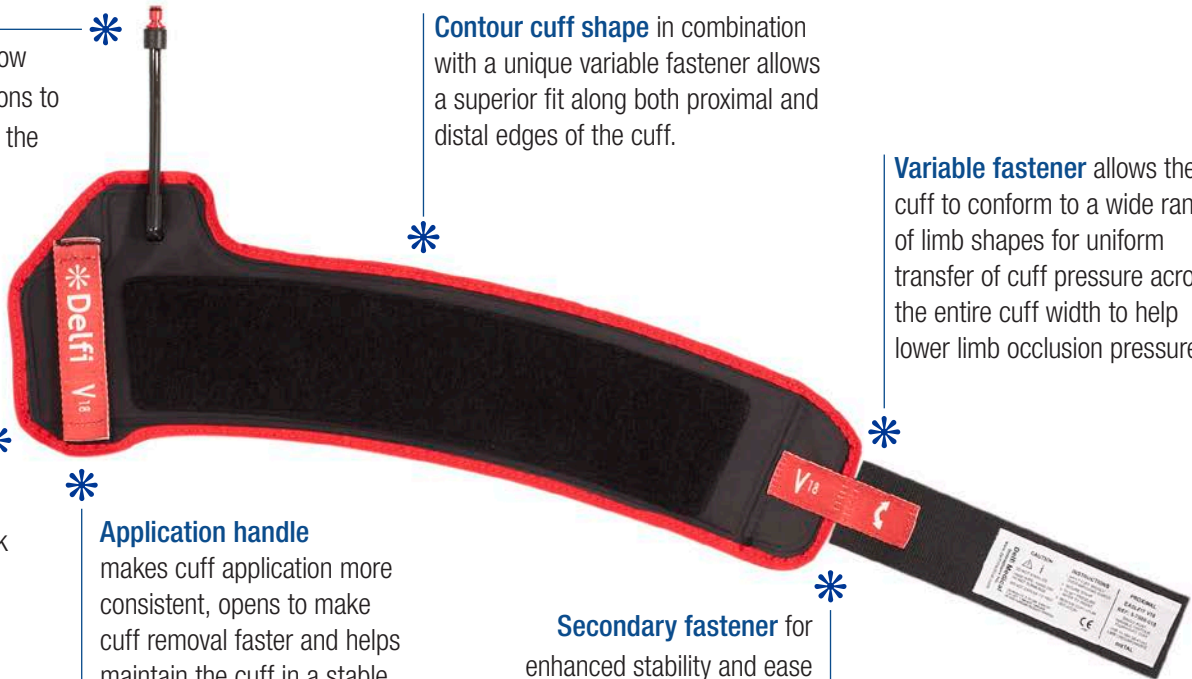

Easi-Fit

Variable Contour and Cylindrical Tourniquet Cuffs

Available in 5 standard lengths to accommodate a variety of limb sizes: 18", 24", 34", 44" and 52" for larger patients.



ORDERING INFORMATION EASI-FIT CUFFS AND MATCHING LIMB PROTECTION SLEEVES (MLPS)

SINGLE PORT		DUAL PORT	MLPS	
9-7500-018	9-7600-018	18" Easi-Fit Variable Contour	9-7920-018	
9-7500-024	9-7600-024	24" Easi-Fit Variable Contour	9-7920-024	
9-7500-034	9-7600-034	34" Easi-Fit Variable Contour	9-7920-034	
9-7500-044	9-7600-044	44" Easi-Fit Variable Cylindrical	9-7920-044	
9-7500-052	9-7600-052	52" Easi-Fit Variable Cylindrical	9-7920-052	

CUFF APPLICATION IS AS EASY AS 1, 2, 3

1 Grasp the application handle and position the ports on the limb at the desired location



2 Slide the securing strap under the handle



3 Tension the cuff for a snug application and secure the strap



CONTACT INFORMATION

P: 800-933-3022 +604-742-0600 **F:** +604-742-3800 **E:** info@delfimedical.com
A: Suite 106-1099 West 8th Avenue, Vancouver, B.C., Canada, V6H 1C3

www.delfimedical.com

U.S. Patent Nos. 5,741,295 6,682,547 7,384,425 7,909,849 and other U.S. and Foreign Patents and Patents Pending

D-4900-010 Rev. 1